


Hands-on: NGS Pipeline

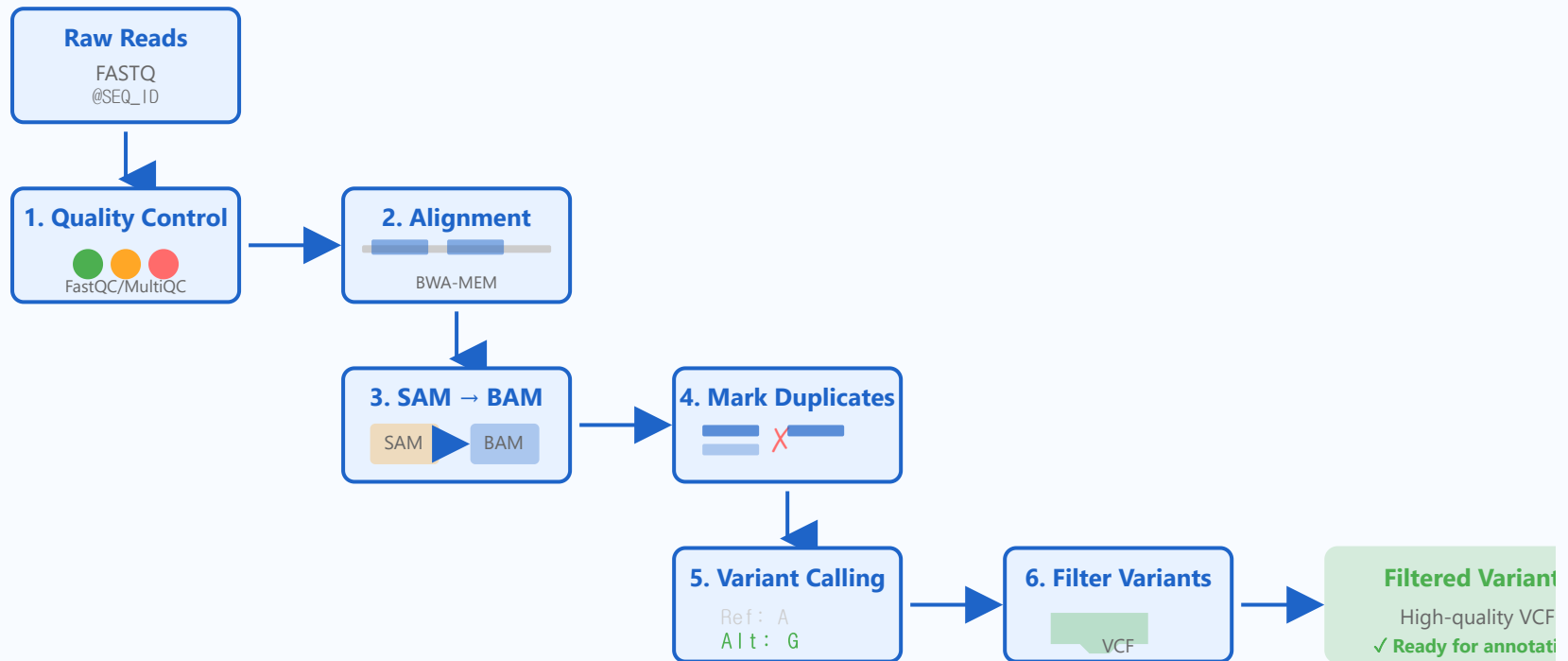

Standard NGS Analysis Pipeline

1. Quality Control

```
fastqc sample_R1.fastq.gz sample_R2.fastq.gz
```

```
multiqc .
```

2. Read Alignment

```
bwa mem -t 8 reference.fa sample_R1.fastq.gz sample_R2.fastq.gz > sample.sam
```

3. Convert SAM to BAM and Sort

```
samtools view -bS sample.sam | samtools sort -o sample.sorted.bam
```

```
samtools index sample.sorted.bam
```

4. Mark Duplicates

```
gatk MarkDuplicates -I sample.sorted.bam -O sample.dedup.bam -M metrics.txt
```

5. Variant Calling

```
gatk HaplotypeCaller -R reference.fa -I sample.dedup.bam -O sample.vcf
```

6. Variant Filtering

```
gatk VariantFiltration -R reference.fa -V sample.vcf -O sample.filtered.vcf
```

Required Software

FastQC, BWA, SAMtools, GATK, Picard

Typical Runtime

4-24 hours depending on coverage and compute resources