


UV-Vis Spectroscopy

$$A = \epsilon bc = -\log_{10}(I/I_0)$$

ϵ : molar absorptivity ($M^{-1}cm^{-1}$) | b : path length (cm) | c : concentration (M)

Spectrophotometer Design

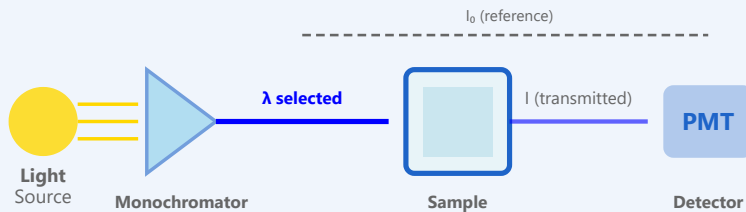

🧬 Chromophores in Biology

Proteins: Trp, Tyr (280 nm)

DNA/RNA: 260 nm

Heme: Soret band (420 nm)

📊 Cuvette Selection

Quartz: UV region

Glass/Plastic: Visible only

Standard: 1 cm path length

🎯 Applications

Protein quantification | DNA/RNA purity | Enzyme kinetics | Drug screening

⚙️ Baseline Corrections

Buffer blank essential | Scatter correction for turbid samples | Temperature control

📈 Linear Range

$A = 0.1$ - 1.0 optimal | Beyond $A=2$: non-linear | Dilute if necessary

💡 Light Sources

Deuterium (UV) | Tungsten-halogen (Visible) | Xenon flash lamps